

The all-new Adapta® Mesa and Adapta® Summit treatment platforms are like no other in function, appearance and value. Take a closer look at the latest technology available, then decide which Adapta platform will be the new foundation you can rely on to treat your patients.

SUMMIT 7

- 7 section treatment platform
- Whisper-quiet hi/lo from 18 to 37 inches controlled by six easy-to-reach cordless, integrated pedal controls around the platform
- Fold away side panels allow easy access to cervical and shoulder region
- Independent positioning of the pivoting head and arm sections
- Gas assist elevates the foot section with the quick-release grab bar. The leg cushion can be raised up to 80 degrees
- Patient positioning system with Adapta® i-Skin standard
- EZ Lift retractable caster system
- PosturFlex option provides flexion/extension and positional traction with the push of a button
- Easily reposition the platform during therapy using the gas-assisted grab bar



ARTICULATIONS




SUMMIT 7

SCHEMATIC

TECHNICAL SPECIFICATIONS

For ultimate versatility in treating all areas of the body. The motorized PosturFlex option "peaks" the center of the table to help accommodate patients who have difficulty lying flat.

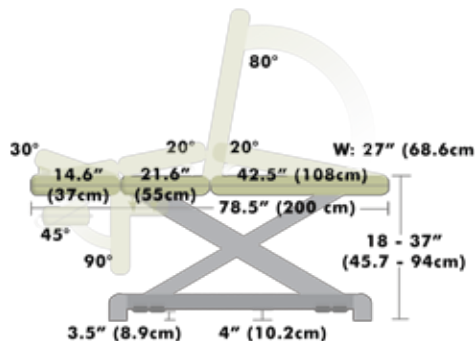

Mains Power: 100-240 V, 50-60 Hz
 Standby Consumption: 7 W
 Input: 230 V / 1 A / 50 Hz ,
 120 V / 1.7 A / 60 Hz
 Weight: 310 lbs. (141 kg)
 Dimensions (L x W x H): 79" x 28" x 18"-37" (200 x 71 x 45.7-94 cm)
 Frame Shipping Dimensions (L x W x H): 60" x 29" x 25" (152 x 74 x 62.5cm)

Cushion Shipping Dimensions (L x W x H): 27" x 6" x 42.5" (69 x 16 x 108cm)
 Lifting Capacity: 450 lbs (204 kg)
 CAL133 Compliant
 Electrical Safety Class: Class II, Type B
 Safety Tests: IEC/UL/EN 60601-1, IEC 60601-1-1, IEC 60601-1-2, CAN/USA C22.2
 No. 601.1-M90 w/ A2



SUMMIT 3

- 3 section treatment platform
- Whisper-quiet hi/lo from 18 to 37 inches controlled by six easy-to-reach cordless, integrated pedal controls around the platform
- Gas-assisted head section pivots up and down independent of the chest section
- Gas assist elevates the foot section with the quick-release grab bar. The leg cushion can be raised up to 80 degrees
- Patient positioning system with Adapta® i-Skin patient surface standard
- EZ Lift retractable caster system
- PosturFlex option provides flexion/extension and positional traction with the push of a button
- Easily reposition the platform during therapy using the gas-assisted grab bar




ARTICULATIONS



SUMMIT 3

A solid platform whose features offer a wide range of treatment options for the busy clinic.

SCHEMATIC

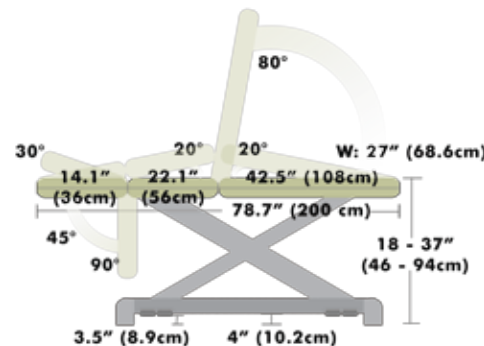

TECHNICAL SPECIFICATIONS

Mains Power: 100-240 V, 50-60 Hz
 Standby Consumption: 7 W
 Input: 230 V / 1 A / 50 Hz,
 120 V / 1.7 A / 60 Hz
 Weight: 310 lbs. (141 kg)
 Dimensions (L x W x H): 79" x 28"
 x 18"-37" (200 x 71 x 45.7-94 cm)
 Frame Shipping Dimensions
 (L x W x H): 60" x 29" x 25"
 (152 x 74 x 62.5cm)

Cushion Shipping Dimensions
 (L x W x H): 27" x 6" x 42.5"
 (69 x 16 x 108cm)
 Lifting Capacity: 450 lbs (204 kg)
 CAL133 Compliant
 Electrical Safety Class: Class II,
 Type B Safety Tests: IEC/UL/EN
 60601-1, IEC 60601-1-1, IEC
 60601-1-2, CAN/USA C22.2
 No. 601.1-M90 w/ A2



The Adapta® i-Skin Patient Surface



- 7 to 9 times more durable than vinyl
- Crack resistant, colorfast, and antimicrobial *
- Non-slip, gently contoured surface helps patients feel secure and comfortable
- Firm enough to keep your elbows and knees from "bottoming out" during manual therapy
- Single piece, solid color molded polyurethane

* Refer to user manual for details

MESA

- 2 section treatment platform
- Whisper-quiet hi/lo from 18 to 37 inches controlled by six easy-to-reach cordless, integrated pedal controls around the platform
- Choose from Naugahyde® upholstery or Adapta® i-Skin patient positioning surface
- Dual face openings allow patients to be oriented on either end to better manage the confines of a small clinic room
- Gas assist elevates one section of the table with the quick-release grab bar. The cushion can be raised to any angle up to 80 degrees
- EZ Lift retractable caster system
- Easily reposition the platform during therapy using the gas-assisted grab bar




ARTICULATIONS



MESA

A functional, durable hi/lo treatment platform for less than what you'd expect.

SCHEMATIC

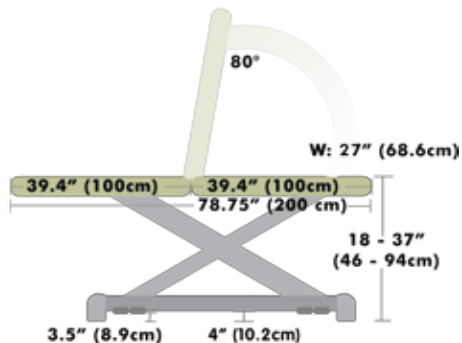

TECHNICAL SPECIFICATIONS

Mains Power: 100-240 V, 50-60 Hz
 Standby Consumption: 7 W
 Input: 230 V / 1 A / 50 Hz ,
 120 V / 1.7 A / 60 Hz
 Weight: 230 lbs. (104 kg)
 Dimensions (L x W x H): 79" x 28"
 x 18"-37" (200 x 71 x 45.7-94 cm)
 Frame Shipping Dimensions
 (L x W x H): 60" x 29" x 25"
 (152 x 74 x 62.5cm)

Cushion Shipping Dimensions
 (L x W x H): 27" x 6" x 42.5"
 (69 x 16 x 108cm)
 Lifting Capacity: 450 lbs (204 kg)
 CAL133 Compliant

Electrical Safety Class: Class II,
 Type B Safety Tests: IEC/UL/EN
 60601-1, IEC 60601-1-1, IEC
 60601-1-2, CAN/USA C22.2
 No. 601.1-M90 w/ A2





Order Details

Part # _____

Color _____

QTY _____

Part # _____

Color _____

QTY _____

Part # _____

Color _____

QTY _____

Part # _____

Color _____

QTY _____

Part # _____

Color _____

QTY _____

Part # _____

Color _____

QTY _____

1.800.592.7329 USA

1.800.242.8329 Fax


1.800.361.6661 Canada

seeadapta.com


a DJO company
DJOglobal.com

Adapta® Summit 7

- 7 section treatment platform
- PosturFlex (optional)
- Adapta® i-Skin patient surface
- 450 lb lift capacity
- Integrated retractable caster system
- Gas-assist moving sections
- 18" – 37" height range adjustable via six easy-to-reach cordless pedal controls

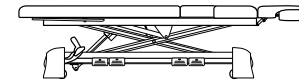

3321 Standard

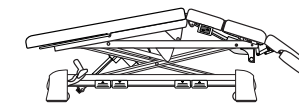


3322 with PosturFlex




Adapta® Summit 3

- 3 section treatment platform
- PosturFlex (optional)
- Adapta® i-Skin patient surface
- 450 lb lift capacity
- Integrated retractable caster system
- Gas-assist moving sections
- 18" – 37" height range adjustable via six easy-to-reach cordless pedal controls

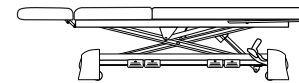

3320 Standard

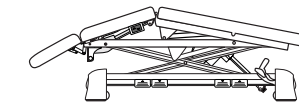


3323 with PosturFlex



Adapta® Mesa

- 2 section treatment platform
- Adapta® i-Skin patient surface
- Economical Naughyde option
- 450 lb lift capacity
- Gas-assist moving sections
- Integrated retractable caster system
- 18" – 37" height range adjustable via six easy-to-reach cordless pedal controls

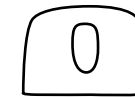


3319 Adapta i-Skin



3317 Naughyde



Naughyde® Color Chart

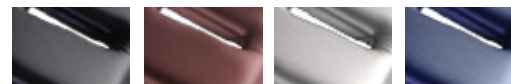

Black Burgundy Dove Imperial Blue

Naughyde® is a registered trademark of Uniroyal Engineered Products, LLC.



Some assembly required. Assembly available upon request for an additional charge.

